

Molecular Cloning

A LABORATORY MANUAL


T. Maniatis Harvard University

E. F. Fritsch Michigan State University

J. Sambrook Cold Spring Harbor Laboratory

W 1 4.27

住友製薬株式会社
研究所図書室

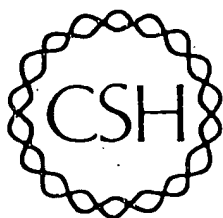

Cold Spring Harbor Laboratory
1982

HYBRIDIZATION OF DNA OR RNA IMMOBILIZED ON FILTERS TO RADIOACTIVE PROBES

There are many methods available to hybridize radioactive probes in solution to DNA or RNA immobilized on nitrocellulose filters. These methods differ in the following aspects:

- the solvent and temperature used (68°C in aqueous solution or 42°C in 50% formamide);
- the volume of solvent and the length of hybridization (large volumes for periods as long as 3 days or minimal volumes for times as short as 4 hours);
- the degree and method of agitation (continuous shaking or stationary);
- the concentration of the labeled probe and its specific activity;
- the use of compounds, such as dextran sulfate, that increase the rate of reassociation of nucleic acids;
- the stringency of washing following the hybridization.

Although the choice depends to a large extent on personal preference, we would like to offer the following guidelines.

1. Hybridization reactions in 50% formamide at 42°C are easier to set up, present less of an evaporation problem, and are less harsh on the filters than is hybridization at 68°C in an aqueous solution. The kinetics of the hybridization reaction in 80% formamide are approximately three to four times slower than in an aqueous solution (Casey and Davidson 1977). Assuming a linear relationship between the rate of hybridization and formamide concentration, the rate in 50% formamide should be two times slower than in an aqueous solution.
2. The smaller the volume of hybridization solvent, the better. The kinetics of nucleic acid reassociation are faster, and the amount of probe needed may be reduced so that the DNA on the filter acts as the driver for the reaction. All these are important parameters when detecting clones of low-abundance mRNAs. However, it is essential that sufficient liquid be present for the filters to remain at all times covered by a film of the hybridization solution.
3. Continual movement of the probe solution across the filter is unnecessary, even for a reaction driven by DNA immobilized on the filter. However, if a large number of filters are hybridized simultaneously, agitation is advisable in order to prevent the filters from adhering to each other.
4. The kinetics of the hybridization reaction are difficult to predict from theoretical considerations, partly because the exact concentration of the immobilized nucleic acid and its availability for hybridization are unknown.

When using probes made by nick translation of double-stranded DNA, the following rule of thumb is useful: Allow the hybridization to proceed for a time sufficient to enable the probe in solution to achieve $1-3 \times C_0t_{1/2}$. In 10 ml of hybridization solution, 1 μ g of a probe of 5-kb complexity will reach $C_0t_{1/2}$ in 2 hours. To determine the time of half-renaturation for any other probe, simply enter the appropriate values into the following equation:

$$\frac{1}{X} \times \frac{Y}{5} \times \frac{Z}{10} \times 2 = \text{number of hours to achieve } C_0t_{1/2}$$

where,

X = the weight of probe added (in μ g)

Y = its complexity (for most probes, complexity is proportional to the length of the probe in kb)

Z = the volume of the reaction (in ml)

After hybridization for $3 \times C_0t_{1/2}$ has been reached, the amount of the probe available for additional hybridization to the filter is negligible. For single-stranded cDNA probes, the hybridization time may be shortened since the lack of a competing DNA strand in solution favors hybridization to DNA bound to the filter.

5. In the presence of dextran sulfate, the rate of association of nucleic acids is accelerated because the nucleic acids are excluded from the volume of the solution occupied by the polymer. Their effective concentration is therefore increased. The rate of association reportedly increases 10-fold in the presence of 10% dextran sulfate (Wahl et al. 1979).

Although dextran sulfate is useful in circumstances where the rate of hybridization is the limiting factor in detecting sequences of interest, it is unnecessary for most purposes. It is also difficult to handle because of its viscosity and sometimes can lead to high backgrounds.

6. In general, the washing conditions should be as stringent as possible; i.e., a combination of temperature and salt concentration should be chosen that is slightly (5°C) below the T_m of the hybrid under study. The temperature and salt conditions can often be determined empirically in preliminary experiments where Southern blots (see pages 382ff) of genomic DNA are hybridized to the probe of interest and then washed under conditions of different stringency.

concentrations of protein are believed to interfere with the annealing of the probe to its target. This quenching of the hybridization signal is particularly noticeable when oligonucleotides or probes less than 100 nucleotides in length are used.

7. In the presence of 10% dextran sulfate or 10% polyethylene glycol, the rate of hybridization is accelerated approximately tenfold (Wahl et al. 1979; Renz and Kurz 1984; Amasino 1986) because nucleic acids are excluded from the volume of the solution occupied by the polymer and their effective concentration is therefore increased. Although dextran sulfate and polyethylene glycol are useful in circumstances where the rate of hybridization is the limiting factor in detecting rare sequences (e.g., northern or genomic Southern blots), they are of no benefit when screening bacterial colonies or bacterial plaques. In addition, they can sometimes lead to high backgrounds, and hybridization solutions containing them are always difficult to handle because of their viscosity. We therefore recommend that dextran sulfate and polyethylene glycol not be used unless the rate of hybridization is very slow, the filter contains very small amounts of DNA, or the amount of radiolabeled probe is limiting.
8. To maximize the rate of annealing of the probe with its target, hybridizations are usually carried out in solutions of high ionic strength ($6 \times$ SSC or $6 \times$ SSPE) at a temperature that is $20-25^{\circ}\text{C}$ below the melting temperature (T_m). Both solutions work equally well when hybridization is carried out in aqueous solvents. However, when formamide is included in the hybridization buffer, $6 \times$ SSPE is preferred because of its greater buffering power.
9. In general, the washing conditions should be as stringent as possible (i.e., a combination of temperature and salt concentration should be chosen that is approximately $12-20^{\circ}\text{C}$ below the calculated T_m of the hybrid under study). The temperature and salt conditions can often be determined empirically in preliminary experiments in which samples of genomic DNA immobilized on filters are hybridized to the probe of interest and then washed under conditions of different stringencies.
10. To minimize background problems, it is best to hybridize for the shortest possible time using the minimum amount of probe. For Southern hybridization of mammalian genomic DNA where each lane of the gel contains $10 \mu\text{g}$ of DNA, $10-20 \text{ ng/ml}$ radiolabeled probe (sp. act. = 10^9 cpm/ μg or greater) should be used and hybridization should be carried out for $12-16$ hours at 68°C in aqueous solution or for 24 hours at 42°C in 50% formamide. For Southern hybridization of fragments of cloned DNA where each band of the restriction digest contains 10 ng of DNA or more, much less probe is required. Typically, hybridization is carried out for $6-8$ hours using $1-2 \text{ ng/ml}$ radiolabeled probe (sp. act. = 10^9 cpm/ μg or greater).
11. Useful facts:
 - a. The T_m of the hybrid formed between the probe and its target may be estimated from the following equation (Bolton and McCarthy 1962):

$$\Rightarrow T_m = 81.5^\circ\text{C} - 16.6(\log_{10}[\text{Na}^+]) + 0.41(\%G + C) - 0.62(\% \text{formamide}) - (600/l)$$

where l = the length of the hybrid in base pairs.

This equation is valid for:

- Concentrations of Na^+ in the range of 0.01 M to 0.4 M. It predicts T_m less accurately in solutions of higher $[\text{Na}^+]$.
- DNAs whose G + C content is in the range of 30% to 75%. Note that the depression of T_m in solutions containing formamide is greater for poly(dA:dT) ($0.75^\circ\text{C}/1\%$ formamide) and less for DNAs rich in poly(dG:dC) ($0.50^\circ\text{C}/1\%$ formamide) (Casey and Davidson 1977).

The equation applies to the "reversible" T_m that is defined by optical measurement of hyperchromicity at OD_{257} . The "irreversible" T_m , which is more important for autoradiographic detection of DNA hybrids, is usually $7\text{--}10^\circ\text{C}$ higher than that predicted by the equation.

Similar equations have been derived for:

- i. RNA probes hybridizing to immobilized RNA (Bodkin and Knudson 1985)

$$T_m = 79.8^\circ\text{C} + 18.5(\log_{10}[\text{Na}^+]) + 0.58(\%G + C) + 11.8(\%G + C)^2 - 0.35(\% \text{formamide}) - (820/l)$$

- ii. DNA:RNA hybrids (Casey and Davidson 1977)

$$T_m = 79.8^\circ\text{C} + 18.5(\log_{10}[\text{Na}^+]) + 0.58(\%G + C) + 11.8(\%G + C)^2 - 0.50(\% \text{formamide}) - (820/l)$$

Comparison of these equations shows that the relative stability of nucleic acid hybrids decreases in the following order: RNA:RNA (most stable), RNA:DNA (less stable), and DNA:DNA (least stable). In aqueous solutions, the T_m of a DNA:DNA hybrid is approximately 10°C lower than that of the equivalent RNA:RNA hybrid. In 80% formamide, the T_m of an RNA:DNA hybrid is approximately 10°C higher than that of the equivalent DNA:DNA hybrid.

- b. The T_m of a double-stranded DNA decreases by $1\text{--}1.5^\circ\text{C}$ with every 1% decrease in homology (Bonner et al. 1978).

The above equations apply only to hybrids greater than 100 nucleotides in length. The behavior of oligonucleotide probes is described in detail in Chapter 11.

For a general discussion of hybridization of nucleic acids bound to solid supports, see Meinkoth and Wahl (1984).